
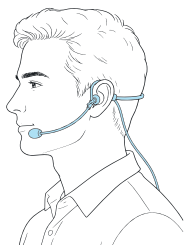


QUICK START


Amplio[™]

Personal Voice Amplification System


- 1 Charge first.** Snap the magnetic cable onto **CHG IN**. The CHG light is red while charging and green when full.
- 2 Wear it around your neck.** Clip the strap to the two bars. Adjust both buckles equally so the amplifier hangs at breastbone level, speaker grille facing your listener.
- 3 Fit the headset.** Band over your ears and around the back of your head, boom on the **left**. Foam windscreen beside the corner of your mouth, one finger-width away.
- 4 Plug in and switch on.** Microphone plug into **MIC IN**. Roll the volume wheel toward ON. The PWR light comes on.
- 5 Set the volume.** Speak with normal effort and raise the volume until listeners hear you clearly. Then back it off slightly: listeners hear more than you do.

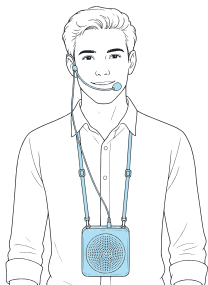

Three rules for clear sound

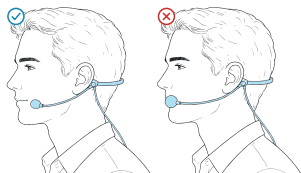

- 1** **Corner of the mouth, never in front.** Breath on the microphone sounds breathy and pops.
- 2** **Never point the speaker at the microphone.** That is what causes the squeal.
- 3** **Use the least volume you need.** If it squeals, turn the volume down, fix the position, then raise it again.

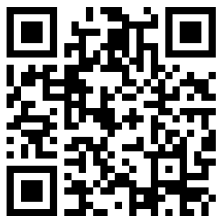

Scan for the full user guide

Setup, wearing, microphones, charging, care, troubleshooting and warranty.

chattervox.store/manuals/amplio

All guides: chattervox.store/manuals

ELT Service Center (800) 286-3481 • service@eltgroup.com

Official ChatterVOX® home: chattervox.store • Nashville, Tennessee, USA

ChatterVOX® is a registered trademark. Keep your sales receipt: the one-year warranty on amplifier and microphone runs from the date of purchase. © 2026 Enhanced Listening Technologies Corp.